

RONIA™ platform

Research Use Only (RUO)

2046-0020 RONIA™ reader

4222-0010 RONIA™ PLGF kit



Complete workflow from whole blood sample to result



Portable reader

- Can run on battery
- No maintenance/service
- No need for laboratory facilities

Quantitative lateral flow test

- PLGF results from finger prick whole blood or serum
- 20 μ L sample volume
- No sample shipping
- Lot specific Calibration for a lot introduced with a RFID


2046-0020 RONIA™ reader



Contents of the case

1. Power adapters
2. USB-C cable
3. Ethernet cable
4. Charger
5. RONIA reader

4222-0010 RONIA™ PIGF kit

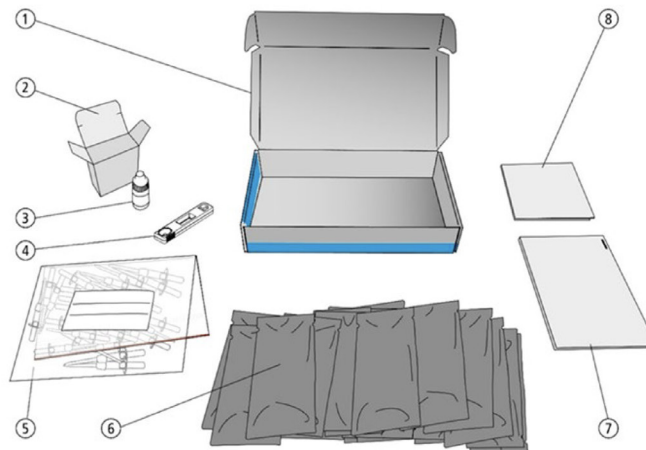

Contents of the kit

1. Kit box
2. Box for Assay Buffer and Calibration Cassette
3. Assay Buffer
4. Calibration Cassette
5. Sample Pipette 20 μ L (28 ea) in a resealable bag
6. Lateral Flow Cassettes (25 ea)
7. Product information sheet
8. Kit insert

Materials required but not supplied with the kit

The following items are required and are available from Wallac Oy or Revvity, Inc., and its distributors.

- Safety-Lancet Extra 18G / Penetration Depth 1.8 mm (prod. no. 85.1017)

In addition, the following materials are required:

- Protective gloves are recommended whenever handling blood or blood products.

- Alcohol swabs for cleaning the fingertips.
- Wipes (used for wiping off the first drop of blood).
- Finger bandages (plasters) may be needed.
- Standard laboratory pipette and tips for 20 μ L if serum is used.

For research use only. Not for use in diagnostic procedures. Please note that product labeling (such as kit insert, product label, and kit box) may be different compared to the company branding. Please contact your local representative for further details.